

Ocean™

WATER SEAL CHEST DRAINS

Traditional Design
you trust...
Innovations
you'll appreciate.

Atrium water seal chest drains set the standard for ease of use, speed of setup, and high performance vacuum control. When it comes to monitoring small changes in patient drainage, fast, accurate fluid assessment has never been easier than with Atrium's improved drainage graduations and highly visible collection chamber designs.

Atrium's fully calibrated water seal provides convenient patient pressure diagnostics and features automatic high negative pressure relief. Best of all, Atrium makes optimal use of a familiar water seal operating system that sets up in seconds and transports easily, safely. Every model is backed with the latest in technology, proven reliability, and unparalleled attention to quality.



The World's Best Selling Chest Drain!

In-Line Connector
with needleless access
sampling port

Internal Knock-Over Nozzles
allow water levels to be recovered in the event of an accidental knock-over


Suction Control Stopcock
provides convenient suction regulation

New Graduated Air Leak Monitor
for assessing patient air leak trends

New Larger Easy-to-Grab Handle and Flexi-Hangers
allow the drain to be hung bedside from a single point

Advanced Float Valve Design
protects against high negative pressure

Large, Easy-to-Read Graphics
provide fast, accurate drainage assessment

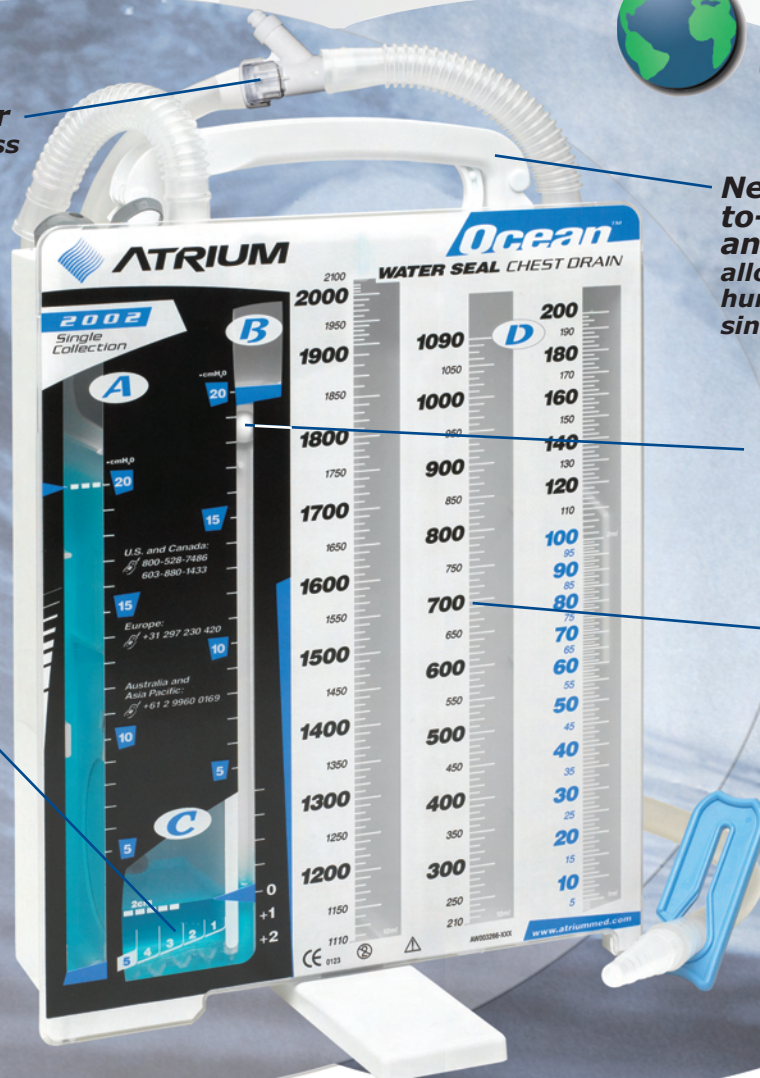

Ocean™

WATER SEAL CHEST DRAINS

Atrium® Ocean™ Water Seal Chest Drains

Ordering Information


2002 Series	Single Collection, 1 Patient Tube	
2002-000	ATS Bag Compatible, Suction Control Stopcock	6
2002-300	Suction Control Stopcock	6
2002-400	Standard System – *Pre-packaged Water	6
2002-050	Pre-connected ATS Bag, Suction Control Stopcock	3
2020 Series	Dual Chamber, 2 Patient Tubes	
2020-000	ATS Bag Compatible, Suction Control Stopcock	6
2020-300	Suction Control Stopcock	6
2012 Series	Infant/Pediatric, 1/4" Tubing, Pediatric Connectors	
2012-320	Suction Control Stopcock	6
2050 Series	ATS Blood Recovery, 1 Patient Tube	
2050-000	ATS Bag Compatible, Suction Control Stopcock	6
2050-300	Suction Control Stopcock	6
2550	In-Line ATS Blood Bag	
	For use with Atrium ATS Compatible Chest Drains	6
2450	Self-Filling ATS Blood Bag	
	For use with Atrium ATS Blood Recovery Drains	10
19911	Drain-Caddy™	
	Chest Drain I.V. Pole Adapter	1



Flexible Options for Post-Op ATS

Whether for routine ATS or for emergency chest trauma applications, Atrium's ATS Blood Recovery models offer unlimited flexibility for prescribing continuous pump reinfusion or post-event blood retrieval with a 2450 Self-Filling ATS blood bag. Models equipped with in-line connectors allow rapid in-line ATS blood bag attachment when required.

Federal (USA) law restricts this device to sale by or on the order of a physician. Please refer to our detailed product catalog for a complete listing of models and pricing. To place an order or to request product information call Atrium toll free at 1-800-528-7486.
*Patented

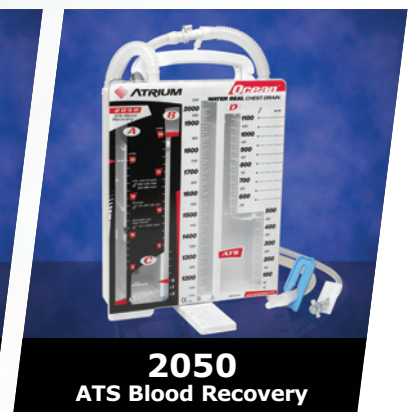

ATRIUM MEDICAL CORPORATION
5 Wentworth Drive
Hudson, New Hampshire 03051 U.S.A.
☎ 603-880-1433 📠 603-880-6718

ATRIUM EUROPE B.V.
Rendementsweg 24 L
3641 SL Mijdrecht, The Netherlands
☎ +31-297-230-420 📠 +31-297-282-653

ATRIUM AUSTRALIA-PACIFIC RIM PTY. LTD.
L 1 Bridgepoint, 3 Brady Street
Mosman NSW 2088 Australia
☎ +61-2-9960-0169 📠 +61-2-8969-2735